




Casero Beds

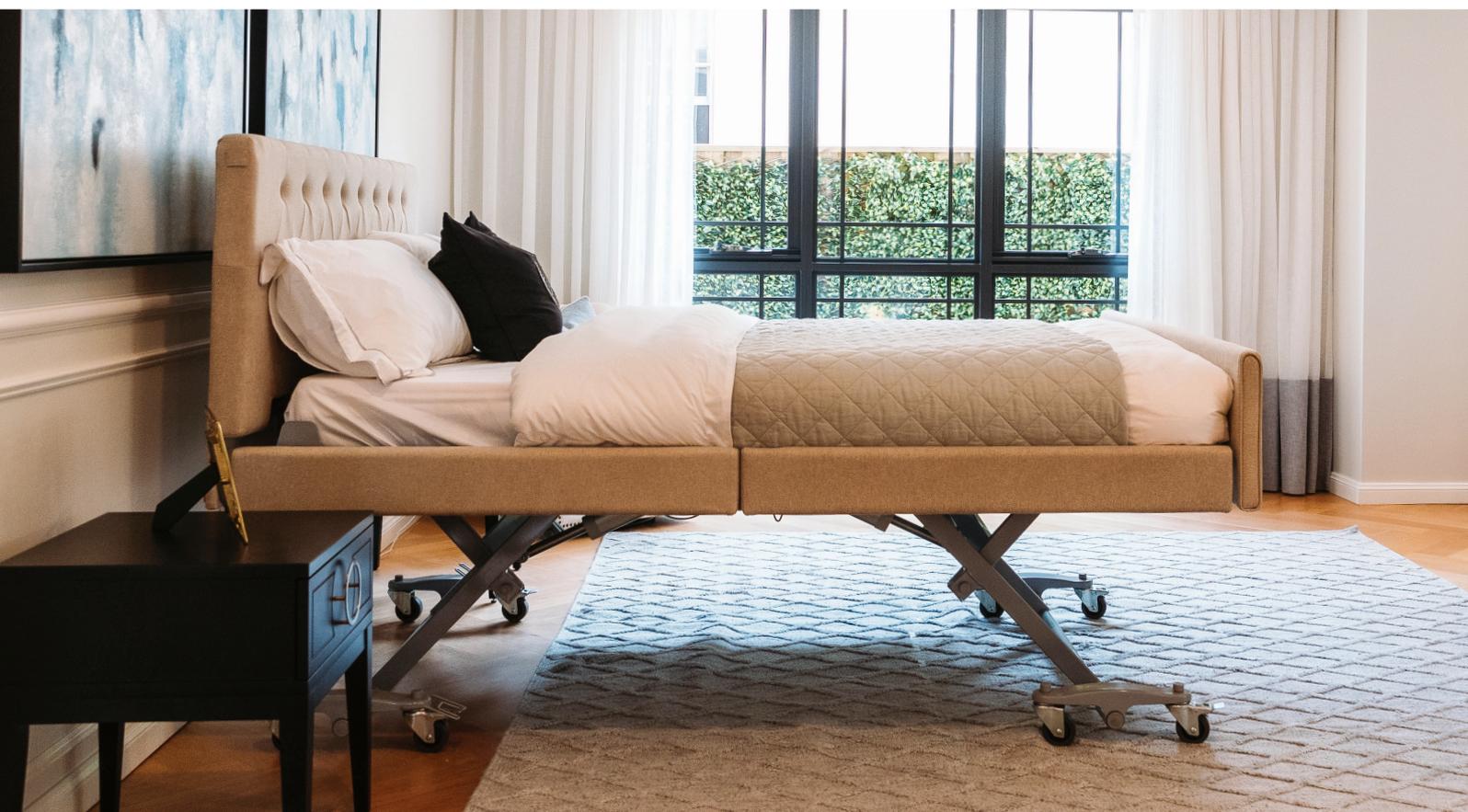

This brochure is for informational purposes only. It may not be reproduced, shared, or stored in any form digital or printed without prior written permission from Movi.

We have taken care to ensure the details in this brochure are accurate at the time of publication; however, specifications, colours, images, and weight limits should be considered approximate and are subject to change without notice.

Movi reserves the right to update or improve products and designs at any time.

Introducing the Movi Casero Bed Range

Comfort. Care. Casero.

The Movi Casero bed range represents a new era in comfort, dignity, and design for home and aged care environments. Thoughtfully rebranded from the trusted Casero range, Movi Casero continues to deliver the same reliability and functionality—now with a refreshed identity that reflects our commitment to innovation and wellbeing.

Designed to blend seamlessly into residential settings, Movi Casero beds offer advanced clinical features without compromising on style. With intuitive controls and advanced functionality, each bed supports both carers and residents in creating a safe, comfortable, and empowering care experience.

Whether for independent living or supported care, Movi Casero is more than a bed—it's a foundation for better living.

Movi Casero Eco Bed



Designed for ease of access in the home environment, the Casero Economy Bed brings reliable functionality and a strong support capacity to everyday care situations. Ideal for the home setting where space, access and simplicity matter, the modular frame delivers the essentials cost-effectively while maintaining the trusted Casero care standard.

Features:

- » Three-piece pull-apart frame ensures the bed can be easily transported and installed through tight doorways, stairwells or restricted access areas without compromising on performance.
- » Stylish head and footboard options to suit aesthetic and infection control preferences.
- » Quality LINAK motors, hand control and electrical components are low maintenance and offer smooth, quiet operation.
- » Electric hi/low functions allow effortless height adjustment for safer transfers and optimal caregiving access.
- » Automatic backrest regression focusses on patient comfort and skin integrity by automatically moving the backrest backward as it elevates, reducing shear and friction, preventing sliding, improving positioning and minimizing sacral pressure.
- » Electric back elevation provides comfortable upper body support, improving breathing, feeding, and overall user comfort.
- » Electric knee break function reduces pressure on the lower body and enhances positional comfort by raising the legs gently.
- » Mesh mattress platform improves breathability for the mattress, assisting microclimate management.
- » Integrated length extension adds up to 200mm to the bed frame, providing extra space and comfort for taller users without compromising stability or functionality.
- » Mattress retainers provide a stable surface for getting into, and out of bed.
- » Individual locking, twin wheel castors improve movability of the bed.
- » Strong build with a safe working load of 200kg, offering dependable support for a wide range of users.
- » Refer to pages 15-16 for range of options and accessories available.

Headboard styles



Colour options



| Specifications | | Single (BD7120) | King Single (BD7130) |
|---|--------------------------|-----------------|----------------------|
|  | Patient Surface Width | 900mm | 1050mm |
| | Patient Surface Length | 2000-2200mm | 2000-2200mm |
|  | Length Extension | 200mm | 200mm |
|  | Overall Width | 1020mm | 1190mm |
| | Overall Length | 2170-2370mm | 2170-2370mm |
|  | Low Height | 280mm | 280mm |
| | High Height | 858mm | 858mm |
|  | Safe Working Load | 200kg | 200kg |
|  | Backrest Range | 0-70° | 0-70° |
|  | Knee Break Range | 0-32° | 0-32° |
|  | Backrest Section Length | 780mm | 780mm |
| | Fixed Section Length | 290mm | 290mm |
| | Upper Leg Section Length | 380mm | 380mm |
| | Lower Leg Section Length | 600mm | 600mm |

Warranty

10 Years Warranty Steel Frame

05 Years Warranty Electrics

01 Years Warranty Hand Control

01 Years Warranty Upholstered surround Including head and foot boards



Scan the QR code
To explore our product configurator and customise your options

Movi Casero Duo Bed



The Movi Casero Duo is a low-height care bed designed to reduce fall risk while supporting safe, independent living. Its expandable width function adapts to changing care needs, providing added space and comfort without needing to replace the bed.

Features:

- » Stylish head and footboard options to suit aesthetic and infection control preferences.
- » Quality LINAK motors, hand control and electrical components are low maintenance and offer smooth, quiet operation.
- » Electric hi/low functions allow effortless height adjustment for safer transfers and optimal caregiving access.
- » Electric back elevation provides comfortable upper body support, improving breathing, feeding, and overall user comfort.
- » Electric knee break function reduces pressure on the lower body and enhances positional comfort by raising the legs gently.
- » Solid mattress platform provides sufficient support for all mattress and sleeping surfaces.
- » Width adjustable bed frame provides greater comfort, safety, and pressure redistribution by adapting to changing user needs, while reducing the need for bed replacement as care requirements evolve.
- » Trendelenburg and reverse Trendelenburg functions improves circulation and supports care procedures, while seat and TV positions provide better postural support.
- » Integrated length extension adds up to 200mm to the bed frame, providing extra space and comfort for taller users without compromising stability or functionality.
- » Central locking, foot-end castors provide safer, simpler and more efficient bed mobility by allowing all wheels to be locked or unlocked with a single foot action.
- » Strong build with a safe working load of 250kg, offering dependable support for a wide range of users.
- » Refer to pages 15-16 for range of options and accessories available.

Headboard styles



Colour options



| Specifications | | Single (BD6000) | King Single (BD6100) |
|---|--------------------------|-----------------|----------------------|
|  | Patient Surface Width | 900mm | 1050mm |
| | Patient Surface Length | 2000mm - 2200mm | 2000mm - 2200mm |
|  | Length Extension | 200mm | 200mm |
|  | Overall Width | 1020mm | 1190mm |
| | Overall Length | 2170mm - 2370mm | 2170mm - 2370mm |
|  | Low Height | 184mm | 184mm |
| | High Height | 760mm | 760mm |
|  | Safe Working Load | 250kg | 250kg |
|  | Backrest Range | 0-65° | 0-65° |
|  | Knee Break Range | 0-32° | 0-32° |
|  | Trendelenburg Range | + - 12° | + - 12° |
|  | Backrest Section Length | 770mm | 770mm |
| | Fixed Section Length | 280mm (Up/Down) | 280mm (Up/Down) |
| | Upper Leg Section Length | 340mm | 340mm |
| | Lower Leg Section Length | 660mm | 660mm |

Warranty

10 Years Warranty Steel Frame

05 Years Warranty Electrics

01 Years Warranty Hand Control

01 Years Warranty Upholstered surround Including head and foot boards



Scan the QR code
To explore our product configurator and customise your options

Movi Casero Mode Bed



The Casero Mode is a floorline bed designed for safe, low-height living, helping reduce falls and injury risks while supporting comfort and dignity at home.

Three-piece modular frame disassembles ensuring the bed can be easily transported and installed through tight doorways, stairwells or restricted access areas without compromising on performance.

Features:

- » Stylish head and footboard options to suit aesthetic and infection control preferences.
- » Quality LINAK motors, hand control and electrical components are low maintenance and offer smooth, quiet operation.
- » Electric hi/low functions allow effortless height adjustment for safer transfers and optimal caregiving access.
- » Automatic backrest regression focusses on patient comfort and skin integrity by automatically moving the backrest backward as it elevates, reducing shear and friction, preventing sliding, improving positioning and minimizing sacral pressure.
- » Electric back elevation provides comfortable upper body support, improving breathing, feeding, and overall user comfort.
- » Electric knee break function reduces pressure on the lower body and enhances positional comfort by raising the legs gently.
- » Solid mattress platform provides support for all mattresses and sleeping surfaces
- » Trendelenburg and reverse Trendelenburg functions improves circulation and supports care procedures, while seat and TV positions provide better postural support.
- » Integrated length extension adds up to 200mm to the bed frame, providing extra space and comfort for taller users without compromising stability or functionality.
- » Central locking, foot-end castors provide safer, simpler and more efficient bed mobility by allowing all wheels to be locked or unlocked with a single foot action.
- » Strong build with a safe working load of 200kg, offering dependable support for a wide range of users.
- » Refer to pages 15-16 for range of options and accessories available.

Headboard styles



Colour options



| Specifications | | Single (BD6200) | King Single (BD6300) |
|---|--------------------------|-------------------------|-------------------------|
|  | Patient Surface Width | 900mm | 1070mm |
| | Patient Surface Length | 2000mm - 2200mm | 2000mm - 2200mm |
|  | Length Extension | 200mm | 200mm |
|  | Overall Width | 1020mm | 1190mm |
| | Overall Length | 2170mm - 2370mm | 2170mm - 2370mm |
|  | Low Height | 110mm | 110mm |
| | High Height | 740mm | 740mm |
|  | Safe Working Load | 200kg | 200kg |
|  | Backrest Range | 0-70° | 0-70° |
|  | Knee Break Range | 0-36° | 0-36° |
|  | Trendelenburg Range | + - 12° | + - 12° |
|  | Backrest Section Length | 880mm | 880mm |
| | Fixed Section Length | 310mm (up) 180mm (down) | 310mm (up) 180mm (down) |
| | Upper Leg Section Length | 330mm | 330mm |
| | Lower Leg Section Length | 530mm - 730mm | 530mm - 730mm |

Warranty

10 Years Warranty Steel Frame

05 Years Warranty Electrics

01 Years Warranty Hand Control

01 Years Warranty Upholstered surround Including head and foot boards



Scan the QR code

To explore our product configurator and customise your options

Movi Casero Vogue Bed



The Movi Casero Vogue is a Floorline bed which significantly reduces the risk of injury from falls while enhancing overall resident safety. Its wide height range and advanced positioning functions support greater range of motion, promoting comfort, mobility, and effective therapeutic positioning for the client.

Features:

- » Stylish head and footboard options to suit aesthetic and infection control preferences.
- » Quality LINAK motors, hand control and electrical components are low maintenance and offer smooth, quiet operation.
- » Electric hi/low functions allow effortless height adjustment for safer transfers and optimal caregiving access.
- » Automatic backrest regression focusses on patient comfort and skin integrity by automatically moving the backrest backward as it elevates, reducing shear and friction, preventing sliding, improving positioning and minimizing sacral pressure.
- » Electric back elevation provides comfortable upper body support, improving breathing, feeding, and overall user comfort.
- » Electric knee break function reduces pressure on the lower body and enhances positional comfort by raising the legs gently.
- » Solid mattress platform provides sufficient support for all mattress and sleeping surfaces.
- » Trendelenburg and reverse Trendelenburg functions improves circulation and supports care procedures, while seat and TV positions provide better postural support.
- » Integrated length extension adds up to 200mm to the bed frame, providing extra space and comfort for taller users without compromising stability or functionality.
- » Central locking, foot-end castors provide safer, simpler and more efficient bed mobility by allowing all wheels to be locked or unlocked with a single foot action.
- » Strong build with a safe working load of 250kg, offering dependable support for a wide range of users.
- » Refer to pages 15-16 for range of options and accessories available.

Headboard styles



Colour options



| Specifications | | Single (BD6400) | King Single (BD6500) |
|---|--------------------------|-------------------------|-------------------------|
|  | Patient Surface Width | 900mm | 1070mm |
| | Patient Surface Length | 2000mm - 2200mm | 2000mm - 2200mm |
|  | Length Extension | 200mm | 200mm |
|  | Overall Width | 1020mm | 1190mm |
| | Overall Length | 2170mm - 2370mm | 2170mm - 2370mm |
|  | Low Height | 110mm | 110mm |
| | High Height | 740mm | 740mm |
|  | Safe Working Load | 200kg | 200kg |
|  | Backrest Range | 0-70° | 0-70° |
|  | Knee Break Range | 0-36° | 0-36° |
|  | Trendelenburg Range | + - 12° | + - 12° |
|  | Backrest Section Length | 880mm | 880mm |
| | Fixed Section Length | 310mm (up) 180mm (down) | 310mm (up) 180mm (down) |
| | Upper Leg Section Length | 330mm | 330mm |
| | Lower Leg Section Length | 530mm - 730mm | 530mm - 730mm |

Warranty

10 Years Warranty Steel Frame

05 Years Warranty Electrics

01 Years Warranty Hand Control

01 Years Warranty Upholstered surround Including head and foot boards



Scan the QR code

To explore our product configurator and customise your options

Movi Casero Plus Bed



The Movi Casero Plus is a bariatric care bed designed for couples or clients needing greater width and weight capacity. Its pull-apart design allows easy home delivery, while the expandable double-to-queen width adapts to changing care needs without replacing the bed.

Features:

- » Two-piece modular frame disassembles ensuring the bed can be easily transported and installed through tight doorways or restricted access areas without compromising on performance.
- » Stylish head and footboard options to suit aesthetic and infection control preferences.
- » Quality motors, hand control and electrical components are low maintenance and offer smooth, quiet operation.
- » Electric hi/low functions allow effortless height adjustment for safer transfers and optimal caregiving access.
- » Electric back elevation provides comfortable upper body support, improving breathing, feeding, and overall user comfort.
- » Electric knee break function reduces pressure on the lower body and enhances positional comfort by raising the legs gently.
- » Solid mattress platform provides sufficient support for all mattress and sleeping surfaces.
- » Trendelenburg and reverse Trendelenburg functions improves circulation and supports care procedures, while seat and TV positions provide better postural support.
- » Integrated length extension adds up to 200mm to the bed frame, providing extra space and comfort for taller users without compromising stability or functionality.
- » Individual locking castors provide safe and simple bed mobility
- » Strong build with a safe working load of 400kg, offering dependable support for a wide range of users.
- » Refer to pages 15-16 for range of options and accessories available.

Headboard styles



Colour options



| Specifications | | Double (BD6600) | Queen (BD6700) |
|---|--------------------------|-----------------|-----------------|
|  | Patient Surface Width | 1350mm | 1500mm |
| | Patient Surface Length | 2000mm - 2200mm | 2000mm - 2200mm |
|  | Length Extension | 200mm | 200mm |
|  | Overall Width | 1490mm | 1630mm |
| | Overall Length | 2170mm - 2370mm | 2170mm - 2370mm |
|  | Low Height | 260mm | 260mm |
| | High Height | 660mm | 660mm |
|  | Safe Working Load | 400kg | 400kg |
|  | Backrest Range | 0-65° | 0-65° |
|  | Knee Break Range | 0-35° | 0-35° |
|  | Trendelenburg Range | + - 12° | + - 12° |
|  | Backrest Section Length | 770mm | 770mm |
| | Fixed Section Length | 280mm (Up/Down) | 280mm (Up/Down) |
| | Upper Leg Section Length | 340mm | 340mm |
| | Lower Leg Section Length | 660mm | 660mm |

Warranty

10 Years Warranty Steel Frame

03 Years Warranty Electrics

01 Years Warranty Hand Control

01 Years Warranty Upholstered surround Including head and foot boards



Scan the QR code
To explore our product configurator and customise your options

Movi Casero Companion Bed - Fixed



The Casero Companion Bed allows partners and caregivers to rest safely and comfortably alongside a Casero bed, supporting connection and peace of mind. Its adjustable height, minimal gap design, and folding frame ensure seamless alignment, easy home delivery and practical everyday use.

Features:

- » Ideal solution for couples, partners, or caregivers who wish to remain close and connected while resting.
- » Provides a safe, stable sleep surface positioned alongside any Casero bed, allowing you to rest comfortably within arm's reach
- » Its thoughtful design minimises the gap between beds, enhancing safety, comfort and interaction.
- » All existing features and functions of the primary Casero bed remain operational when paired as a Companion Bed set up.
- » The Companion Bed is available as a manual or electric adjustable frame and can be configured on either side of the primary bed for maximum flexibility.
- » The durable steel frame offers multiple height positions, ensuring alignment with the primary bed for a companion setup.
- » Bed frame folds in half to assist deliveries or for ease of storage when not in use.
- » Compatible with Casero Eco, Duo, Vogue and Mode beds, in Single and King Single sizes.
- » Refer page 15-16 for range of options and accessories available

Headboard styles



Colour options



| Specifications | | Single (BD6830) | King Single (BD6840) |
|----------------|------------------------|-----------------|----------------------|
| | Patient Surface Width | 880mm | 1050mm |
| | Patient Surface Length | 1980mm | 1980mm |
| | Length Extension | 200mm | 200mm |
| | Overall Width | 1020mm | 1190mm |
| | Overall Length | 2170mm - 2370mm | 2170mm - 2370mm |
| | Low Height | 140mm | 140mm |
| | High Height | 300mm | 300mm |
| | Safe Working Load | 200kg | 200kg |

Movi Casero Companion Bed - Electric



The Casero Electric Companion Bed enables seamless alignment with all Casero beds—Eco, Duo, Vogue and Mode—allowing couples to rest together safely and comfortably. With independent electric bed adjustments and precise positioning, it supports true side-by-side alignment while maintaining individual care and independence.

Features:

- » Electric height adjustment allows precise alignment with Casero Eco, Duo, Vogue and Mode beds.
- » Supports safe, side-by-side sleeping for couples while maintaining individual care needs.
- » It's considered design minimises the gap between beds, enhancing safety, comfort and interaction.
- » The Companion Bed is available as a manual or electric adjustable frame and can be configured on either side of the primary bed for maximum flexibility.
- » Durable, stable construction suited to everyday home and care environments.
- » Refer page 15-16 for range of options and accessories available

Warranty

| | |
|----|--|
| 10 | Years Warranty Steel Frame |
| 05 | Years Warranty Electrics |
| 01 | Years Warranty Hand Control |
| 01 | Years Warranty Upholstered surround Including head and foot boards |

| Compatible Electric Companion Beds | |
|------------------------------------|-------------------------------------|
| Eco Bed | Single - two beds side by side |
| Eco Bed | King Single - two beds side by side |
| Duo Bed | Single - two beds side by side |
| Duo Bed | King Single - two beds side by side |
| Mode Bed | Single - two bed side by side |
| Mode Bed | King Single - two beds side by side |
| Vogue Bed | Single - two beds side by side |
| Vogue Bed | King Single - two beds side by side |



Scan the QR code
To explore our product configurator and customise your options

Movi Casero Bed

Accessories & Options




| Movi Casero Beds | Single | King Single |
|------------------------------------|-----------------|-----------------|
| Movi Casero Eco Bed | BD7120 | BD7130 |
| Movi Casero Duo Bed | BD6000 | BD6100 |
| Movi Casero Mode Bed | BD6200 | BD6300 |
| Movi Casero Vogue Bed | BD6400 | BD6500 |
| Movi Casero Fixed Companion Bed | BD6830 | BD6840 |
| Movi Casero Electric Companion Bed | BD6830 + BZ0660 | BD6840 + BZ0660 |
| Movi Casero Beds | Double | Queen |
| Movi Casero Plus Bed | BD6600 | BD6700 |


Headboard styles




Footboard styles

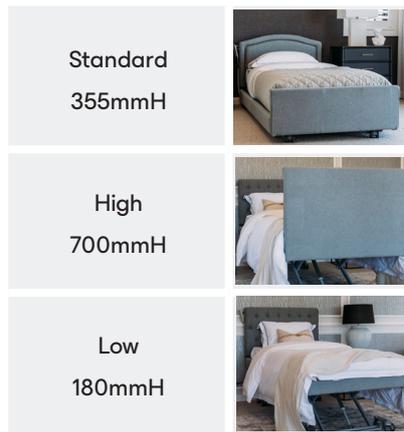


Colour options

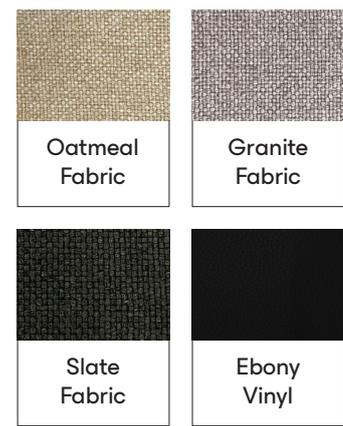



Padded Panels

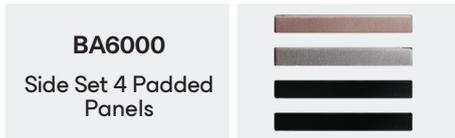


Side Options

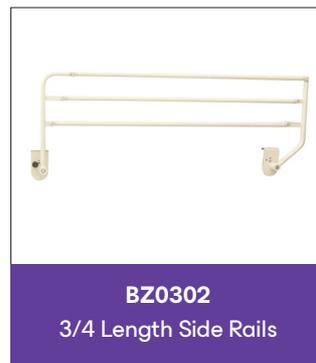

Accessories

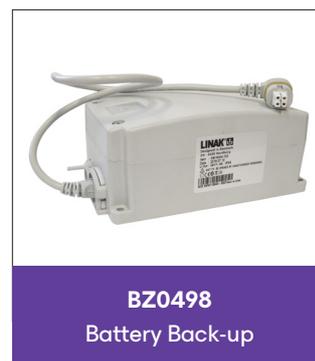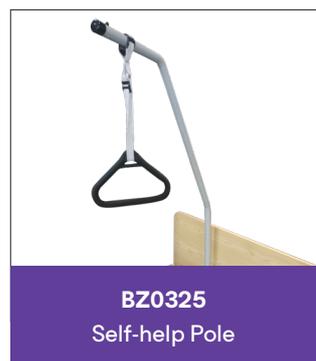

Movi Casero Luxury Comfort Mattresses

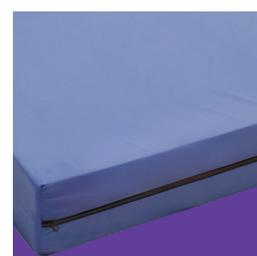

- » Three-layer foam construction with a memory-foam top comfort layer delivers pressure relief and contours to body shape
- » Firm core/side support structure designed to assist with transfers and maintain edge stability
- » Two cover material options: bamboo-quilted home-care style cover or a medical-grade water resistance cover, offering flexibility depending on care environment.
- » Specifically engineered for low-level pressure-care and medium-level comfort — ideal for users requiring support and durability without overly plush feel.

Specifications

| | Size | Code | Cover | Feel | Height |
|-------------------------------|-------------|--------|---------------|--------|--------|
| Movi Casero Luxury Comfort | Single | MF4010 | Bamboo | Medium | 150mm |
| | King Single | MF4020 | Bamboo | Medium | 150mm |
| | Double | MF4030 | Bamboo | Medium | 150mm |
| | Queen | MF4040 | Bamboo | Medium | 150mm |
| | Single | MF4110 | Pressure Care | Medium | 150mm |
| | King Single | MF4120 | Pressure Care | Medium | 150mm |
| | Double | MF4130 | Pressure Care | Medium | 150mm |
| | Queen | MF4140 | Pressure Care | Medium | 150mm |

Movi Casero Luxury Comfort Plus Mattresses



- » Four-layer foam construction features a memory-foam top for contouring and pressure redistribution.
- » Firm comfort profile with a supportive core and reinforced edge/sides to aid safe transfers and maintain mattress stability
- » Two cover material options: bamboo-quilted home-care style cover or a medical-grade water resistance cover, offering flexibility depending on care environment.
- » Designed for low-level pressure care needs while delivering durable, long-term support in home or facility settings.

| Specifications | | | | | |
|------------------------------------|-------------|--------|---------------|------|--------|
| | Size | Code | Cover | Feel | Height |
| Movi Casero Luxury Comfort Plus | Single | MF4210 | Bamboo | Firm | 200mm |
| | King Single | MF4220 | Bamboo | Firm | 200mm |
| | Double | MF4230 | Bamboo | Firm | 200mm |
| | Queen | MF4240 | Bamboo | Firm | 200mm |
| | Single | MF4310 | Pressure Care | Firm | 200mm |
| | King Single | MF4320 | Pressure Care | Firm | 200mm |
| | Double | MF4330 | Pressure Care | Firm | 200mm |
| | Queen | MF4340 | Pressure Care | Firm | 200mm |



LIVE BETTER